



StoneWrapper®

The high performer of stone extraction baskets...

"StoneWrapper®"



Perfect for removing concretions of different sizes from hard-to-reach areas of the renal pelvis system comfortably, safely and quickly.

- ♦ Perfect hybrid of grabber and stone extraction basket
- ♦ Elastic, fully articulated three-armed nitinol basket
- Easy insertion into the lower calyx and quick removal of all concretions
- Available in two sizes and with two opening options

Effective: work efficiently, save time

The innovative StoneWrapper solves the major challenges involved with flexible ureterorenoscopy: fast and safe atraumatic removal of marginal stones from the renal pelvic system.

The combination of

- a specially designed three-arm hybrid head (combination of basket and grabber),
- a new Guide wire technology with high shaft stiffness and hyperflexible distal end,
- including a Tuohy-Borst adapter and an insertion aid for easy insertion into the endoscope's working channel and as splash guard,
- a wire design designed for durability and functionality,
- an ergonomic push handle for fatigue-free working and
- excellent irrigation for best visibility enables the reliable and fast removal of even the smallest concretions from any position

Size (FR/CH)	Length (cm)	Opening Width Ø (mm)	Order Code
1.5	120	8	48012553
1.5	120	11	48012554
2.0	120	8	48012555
2.0	120	11	48012556


Distributed by UROTECH.



STONE BASKETS

High elastic Nitinol material, with detachable handle acc. to Lahme

"TIPLESS"

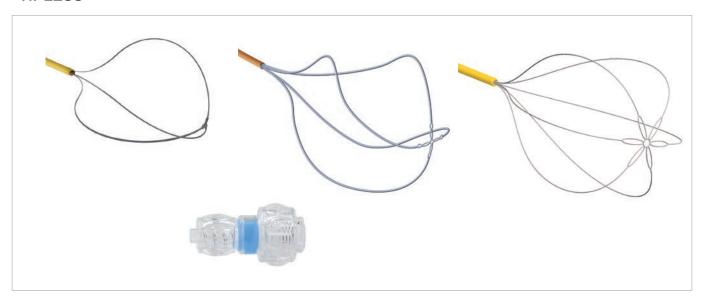

- ♦ Delicate 3-, 4- or 6-wire-baskets without tip for the flexible uretero-renoscopy
- ♦ For grasping of smallest stone concrements also on tissue edge
- With Tuohy Borst Adapter to protect the basket tip

Size (FR/CH)	Wires	Working Length (cm)	Basket Length (mm)	Basket Ø (mm)	Basket Shape	Order Code
1.2	3	120	12	8	straight	STN-T123-120
1.8	4	120	15	11	straight	STN-T184-120
1.8	4	120	15	11	helical	STN-HT184-120
1.8	6	120	13	9	straight	STN-T186-120
2.2	4	120	18	11	straight	STN-T224-120

Handle acc. to Lahme for all stone baskets

- Dis- and reconnectable handle: endoscope can be removed while basket is situated.
 Handle and basket can get attached afterwards
- Fixation of the stone in the basket due to spring mechanism
- ♦ Basket opens and closes within image area: no adjustment of endoscope necessary



"HELICAL"




Size (FR/CH)	Wires	Working Length (cm)	Basket Length (mm)	Basket Ø (mm)	Tip length (mm)	Order Code
2.5	4	90	30	12	3	STN-L254-090
2.5	4	120	30	12	3	STN-L254-120
3.0	4	90	31	12	3	STN-L304-090
3.0	4	120	31	12	3	STN-L304-120
4.0	4	90	37	15	3	STN-L404-090

"FILIFORM"



Size (FR/CH)	Wires	Working Length (cm)	Basket Length (mm)	Basket Ø (mm)	Tip length (mm)	Order Code
3.0	4	90	30	12	40	STN-F304-090

"STRAIGHT"

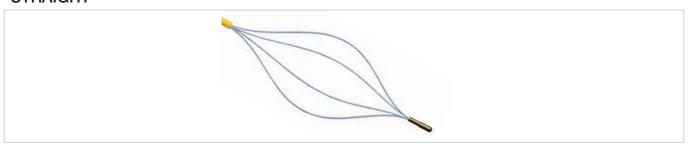

Size (FR/CH)	Wires	Working Length (cm)	Basket Length (mm)	Basket Ø (mm)	Tip length (mm)	Order Code
2.5	4	120	30	12	3	STN-G254-120
3.0	4	120	31	12	3	STN-G304-120